







MODEL 9000A

SMARTSINE[™] CALIBRATION DRIVER

- Acceleration, velocity, displacement, and acoustic readout units
- Battery-powered unit is self-contained solution
- Compact and integrated sinusoidal source, ICP inputs, and meter system
- Offers sinusoidal calibration closed loop control over a 0.5 Hz to 20 kHz frequency range
- Portable by design with heavy duty case

TYPICAL APPLICATIONS

- Accelerometer or velocity sensor calibration
- Large sensor verification (Geophones, Vibration Switches, Seismoprobes, etc.)
- Acoustic calibration: drive acoustic couplers for "back-toback" calibration of array microphones
- Testing of Electronic devices with embedded vibration measurement capability

DRIVE SHAKERS & CREATE CALIBRATION CERTIFICATES

The SmartSine[™] Calibration Driver Model 9000A is a batterypowered sinusoidal signal source used to drive shaker systems to predefined vibration levels and frequencies to calibrate individual sensors, vibration switches and data collectors. The 9000A can also be used to drive an acoustic coupler for calibration of array microphones.

Model 9000A offers two sensor input channels and one output channel. The output channel provides signal input to drive shaker systems and other excitation sources over a wide amplitude and adjustable frequency range. Built-in CALROUTE functionality adds semi-automation capability with pre-programmed frequency and amplitude points using Microsoft Excel[®].

The 9000A displays test sensor sensitivity on the readout screen in real-time by comparing output against a known reference sensor. Both channels feature ICP[®] (IEPE) inputs for common piezoelectric accelerometer and microphone signal conditioning. In addition, users can save up to 500 calibration records directly to the unit's internal memory. Records are transferable via the unit's USB port to a flash drive and imported as an Excel spreadsheet, allowing the creation of ISO 17025-compliant customizable calibration certificates on a computer.

SPECIFICATIONS					
Performance					
Test Sensor In					
Sensor Type	Voltage or ICP® (IEPE) [1]				
Input Voltage (Max)	5 V pk ^[2]				
Reference Sensor In	1				
Sensor Type	ICP® (IEPE) [1] ONLY				
Input Voltage (Max)	1 V pk [2]				
Reference Sensitivity	Single point or curve (up to 60 points)				
Bias Fault Indication (ICP® Sensors)	Yes				
Source Out	Sine wave form of 1 V RMS amplitude Max ^[3]				
Monitor Reference Out	Buffered reference output				
Electrical Connector (Input/Output)	BNC Jack				
Operating Range [4]	0.5 Hz-20 kHz (30-1200k CPM)				
Units of Readout					
Acceleration (pk and RMS)	[g], [m/s ²]				
Velocity (pk and RMS)	[in/s], [mm/s]				
Displacement (pk to pk)	[mils], [µm]				
Sound Pressure Level [5]	[dB]				
Frequency	Hz, (CPM)				
Test Sensor Sensitivity	mV/EU [6]				
Shaker Displacement Limit Setting (Max) ^[7]	10 in pk-pk (254 mm pk-pk)				
Storage	Up to 500 calibration records				
Points Per Record	30 calibration data points				
Sensor Information	Model number, serial number, mounting orientation (x, y, z), and user notes				
Export File Format	Export calibration records to USB flash drive (FAT 32) in CSV (comma-separated values) format				
AC Power (for recharging battery)	110 V–240 V, 50 Hz–60 Hz				
Input Power Rating from Charger	18 VDC, 1 A				
Internal Battery (Li-ion)	11.1 VDC, 2200 mAh				
Operating Battery Life [8]	16 hours				
Physical					
Dimensions (H x W x D)	4.7 x 9.8 x 11.8 in (12 x 25 x 30 cm)				
Weight	7.8 lb (3.56 kg)				
Operating Temperature	32 °F to 122 °F (0 °C to 50 °C)				

[1] 5 mA constant current excitation to ICP® (IEPE) reference sensor

[2] Recommended SUT voltage > 20 mV pk; REF voltage > 4 mV pk

[3] 50 Ω output impedance

[4] Or limited to frequency range of reference sensitivity curve input
[5] Reference 20 μPa

[6] EU: [g], [m/s²], [in/s], [mm/s], [mils], [µm], or [Pa] for acoustics

[7] Actual Max. amplitude may be limited by exciter specifications and amplifier gain

[8] As shipped from factory in new condition

Supplied Accessories [8]

Accessory Pouch Containing:


USB Flash Drive with Calibration Report Generation Worksheet

Universal Power Supply (Power Charger) with interchangeable plug adaptors

Ref Sensitivity	Vibration					
(mV/g)	g pk	in/s pk	mils pk-pk	m/s² pk	mm/s pk	µm pk-pk
10	100	61.45	195.59	981	1561	4968
100	10	6.14	19.56	98	156	497

Ref Sensitivity	Sound Pressure Level
(mV/Pa)	dB [5]
12.5	130
50	115

Maximum Testing Amplitudes

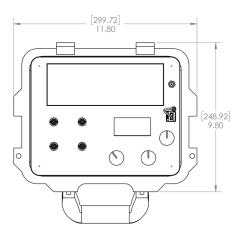

Model 9000A Technical Drawing Dimensions in inches [mm]



10310 Aerohub Boulevard, Cincinnati, OH 45215 USA

modalshop.com | info@modalshop.com | 800 860 4867 | +1 513 351 9919

© 2021 PCB Piezotronics - all rights reserved. PCB Piezotronics is a wholly-owned subsidiary of Amphenol Corporation. Endevco is an assumed name of PCB Piezotronics of North Carolina, Inc., which is a wholly-owned subsidiary of PCB Piezotronics, Inc. Accumetrics, Inc. and The Modal Shop, Inc. are wholly-owned subsidiaries of PCB Piezotronics, Inc. The Modal Shop, Inc. are wholly-owned subsidiaries of PCB Piezotronics, Inc. Ind Sensors and Larson Davis are Divisions of PCB Piezotronics, Inc. Except for any third party marks for which attribution is provided herein, the company names and product names used in this document may be the registered trademarks or unregistered trademarks of PCB Piezotronics, Inc., PCB Piezotronics of North Carolina, Inc. (d/b/a Endevco), The Modal Shop, Inc. or Accumetrics, Inc. Detailed trademark ownership information is available at www.pcb.com/trademarkownership.